






興田國際集團
HTE INTERNATIONAL GROUP



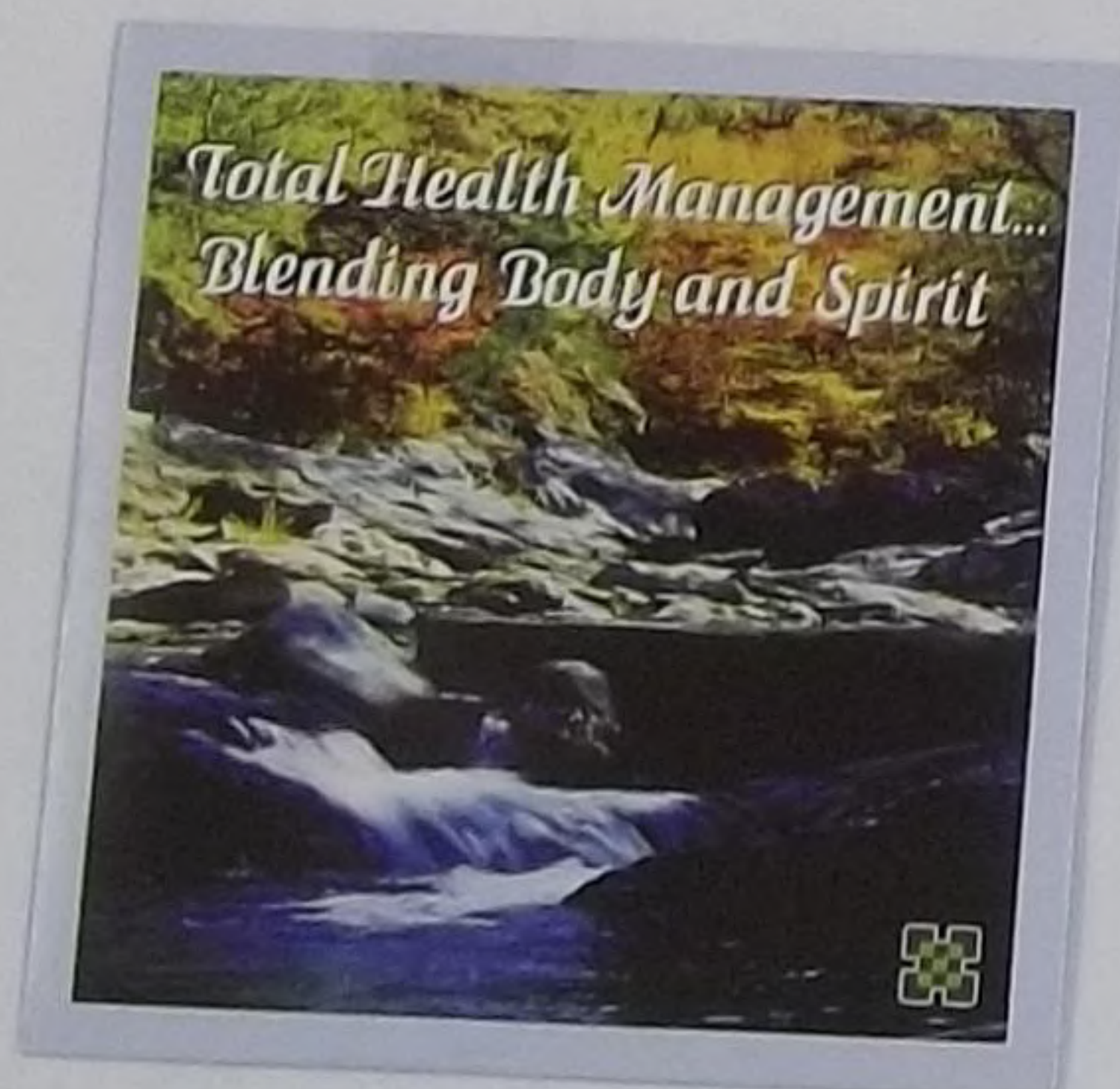



TIMER

OFF

II

30=
min..



=10

20

★ TURN THE DIAL MORE
THAN A QUARTER AND
BACK TO YOUR DESIRED
TIME WHEN SETTING IT
WITHIN FIVE MINUTES



SUN ANCON[®]

SUN ANCON[®] SDM-888

AEROBIC EXERCISER
RATED VOLTAGE: 120V 60Hz
POWER CONSUMPTION: 34W
RATED TIME: 30 MINUTES
LOT NO.: DEG2000467





 **HTE INTERNATIONAL GROUP**
TECHNICAL COLLABORATION SKYLITE CORPORATION, JAPAN

SUN ANCON[®] (SUN HARMONY[®])
THE CHI MACHINE

SDM-888
USA PAT. NO: 5107822
U.K. PAT. NO: 2023294
GERMANY NO: G9208487.7
CHINA NO: ZL201030168135.1
TAIWAN NO: 58736
THAILAND NO: 4056
Specifications: 120V/60Hz/34W
HTE INTERNATIONAL GROUP
Technical collaboration: Skylite corporation Japan

爽安康[®] SDM-888
美國專利文號: 5107822
英國專利文號: 2023294
德國專利文號: G9208487.7
中國專利文號: ZL201030168135.1
台灣專利文號: 58736
泰國專利文號: 4056
規格: 120V/60Hz/34W
興田國際集團
(由興田國際集團與日本史丹利株式會社合作生產)

SERIAL NO.
製造號碼: DEG2000467

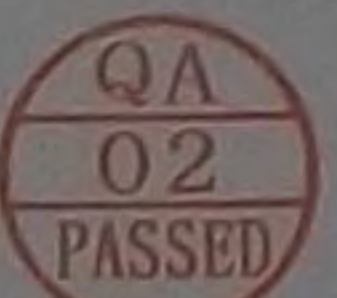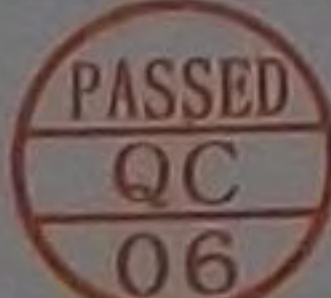生管 QA		品管 QC	
日期 DATE	2011-12-15	日期 DATE	2011-12-15



DEG2000467



SUN ANCON®

SDM-888

AEROBIC EXERCISER

RATED VOLTAGE: 120V 60Hz

POWER CONSUMPTION: 34W

RATED TIME: 30 MINUTES

LOT NO.: DEG2000467



HTE INTERNATIONAL GROUP

TECHNICAL COLLABORATION SKYLITE CORPORATION JAPAN

SUN ANCON® (SUN HARMONY®)

THE CHI MACHINE

SDM-888

USA PAT. NO:

5107822

U.K. PAT. NO:

2023294

GERMANY NO:

G9208487.7

CHINA NO:

ZL201030168135.1

TAIWAN NO:

58736

THAILAND NO:

4056

Specifications: 120V 60Hz/34W

HTE INTERNATIONAL GROUP

Technical collaboration: Skylite corporation Japan

奕安康®

SDM-888

美國專利文號:

5107822

英國專利文號:

2023294

德國專利文號:

G9208487.7

中國專利文號:

ZL201030168135.1

台灣專利文號:

58736

泰國專利文號:

4056

規格: 120V 60Hz/34W

興田國際集團

(由興田國際集團與日本史丹尼特株式會社合作生產)

SERIAL NO.

製造號碼: DEG2000467

生管
QA



品管
QC



日期
DATE

2011-12-15

日期
DATE

2011-12-15



DEG2000467



KEEP YOUR HAND AWAY FROM HERE

PADA SAAT MENGGUNAKAN ALAT KE SEHATAN, JAUHKAN TANGAN ANDA DARI MESIN

ชั้น อัน คั้น ขอบ สุขภาพที่ดีให้ยกถาด

使用時 手請遠離此處



●KEEP YOUR HANDS AWAY FROM THE MOVING ANKLE REST.



●DO NOT STAND ON THE ANKLE REST.







PLACE YOUR FEET BELOW THE
ANKLES ON THIS LEG REST.

LETAKKAN KEDUA KAKI ANDA
DI ATAS TUMPUAN KAKI DEN-
GAN PERLAHAN

โปรดวางข้อเท้าลง
บนถ้ำอย่างช้าๆ

請將雙腳輕放於腳臺上




PLACE YOUR FEET BELOW THE
ANKLES ON THIS LEG REST.

LETAKKAN KEDUA KAKI ANDA
DI ATAS TUMPUAN KAKI DEN-
GAN PERLAHAN

ไปวางข้อเท้าลง
บนเตียงอย่างช้าๆ

請將雙腳輕放於腳臺上

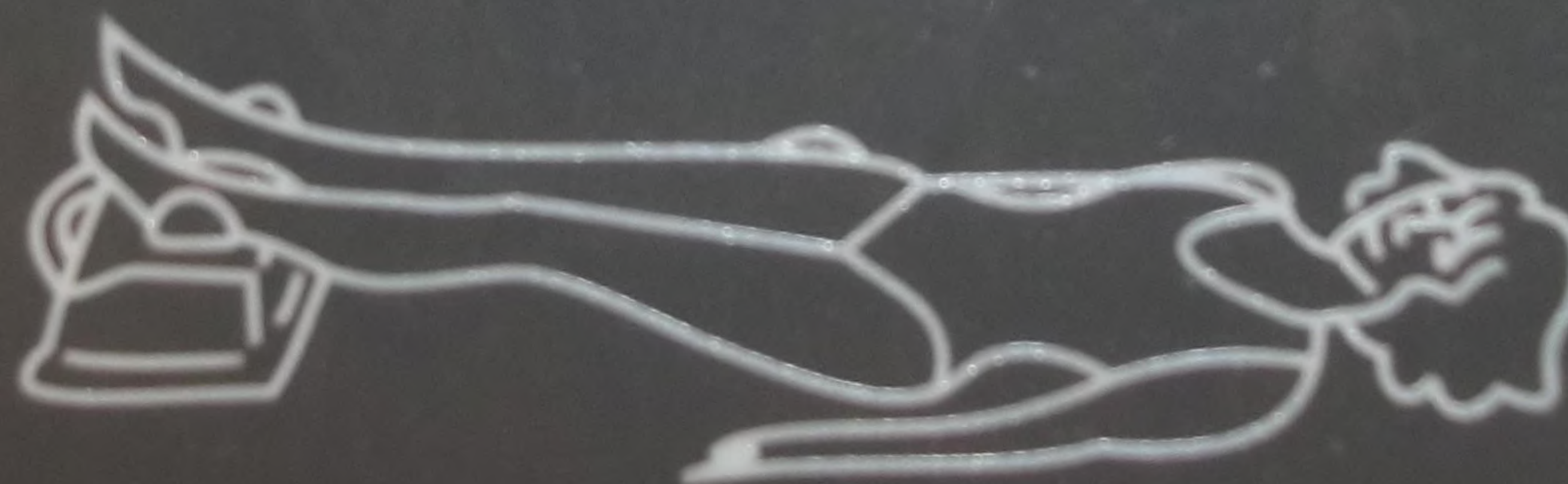







THE CHI MACHINE®

使用指南

- USER'S GUIDE
- PEMAKAIAN
- 사용설명서
- แนะนำการใช้
- MANUEL DE L'UTILISATEUR
- GUIA DEL USUARIO
- BENUTZERHANDBUCH



CAUTIONS:

Do not use the machine within 60 minutes of eating a meal or drinking alcohol. Consult your doctor before using this machine if you suffer from any of the following:

- (1) Any malignancies (such as cancer)
- (2) Serious cardiac ailments (such as a recent heart attack)
- (3) Fractured bones or have had an operation within the past three months, or are under long-term care.
- (4) Serious infection or bleeding injuries or epilepsy.

Start out by limiting the duration of a session to two to three minutes. Increase the duration very gradually if you suffer from minor heart conditions.

Do not use the machine when close to full term pregnancy or if you are known to be prone to spontaneous miscarriage.

If there is any pain, discontinue use until the pain subsides.

WARNING:

When your ankles are correctly placed on the footrest they should feel fairly comfortable both before and after the machine is turned on. However, if you are very sensitive and your ankles still feel uncomfortable even after adjusting their position on the footrest, we suggest placing a small towel under them. Using a towel somewhat reduces the beneficial massage in the ankle area. Unless you are very sensitive, it is better to use the machine without a towel.

WARNING: Do not use the machine in any of the following manners:

1. In an inclined or otherwise unstable position
2. Connected to an electrical receptacle shared by other appliances
3. On a blanket or fuzzy surface (blocking the air vents can cause overheating and malfunction).
4. At temperatures below 0°C degrees (32°F) or above 40°C degrees (104°F), or in areas of high humidity or moisture.
5. If the power supply cord is damaged, send the unit back to HTE and have it replaced by our repair department to prevent an electrical hazard.
6. Pull out the plug while not in use.

USER MAINTENANCE & STORAGE:

CLEANING: Wipe with soft, damp cloth and mild detergent. Be sure to rinse and dry it with a soft cloth.

SERVICING: This unit contains NO USER-SERVICABLE PARTS. Return to the manufacturer or authorized representative for any and all servicing.

STORAGE: When not in use, store the machine in a dry place out of direct sunlight. Take care not to damage the power supply cord or the cord to the timer --avoid bending, twisting or stressing cords at their connectors.

REPAIR SERVICES FREE OF CHARGE

1. Product can be repaired free of charge if it is within the warranty period and if user treats the device with proper care.
2. Damage caused by unauthorized disassembly or operation in contravention to these directions is not covered under the free repair warranty.
3. Free repair is NOT provided even during the warranty period if:
 - (a) The machine has been damaged due to dropping
 - (b) The machine was damaged by acts of God such as earthquake, public disturbances, fire, etc.
4. Please consult the company for other situations, if any.

OWNER WARRANTY REGISTRATION CARD

The Chi Machine® (Sun Ancon Chi Machine)

OWNER'S NAME MR..

MS..

STREET
ADDRESS

CITY

STATE

ZIP

PHONE NO.()

POST CODE

DATE OF PURCHASE

INVOICE NO.

MACHINE SERIAL NO.

SPONSOR'S

SPONSOR'S

NAME

ID NO.

PHONE NO.()

IMPORTANT:

● PLEASE COMPLETE AND MAIL THIS FORM RIGHT AWAY. IT WILL HELP ENSURE THAT H.T.E HAS ACCURATELY REGISTERED YOUR PURCHASE FOR WARRANTY.

● IT WILL ALSO HELP US CONTACT YOU IN THE UNLIKELY EVENT OF A PRODUCT RECALL.

History

The Chi Machine® Therapeutic Massager helps virtually anyone locally relax muscles with little or no conscious effort. Just recline, face-up on the floor with both feet placed in the specially designed footrest. The Chi Machine® rocks the entire body from side to side in a gentle undulating motion.


The Chi Machine® Therapeutic Massager is the culmination of 38 years of research and development by Dr. Shizuo Inoue, a medical visionary and the developer of the world renowned Chi Machine. The patented precision design of this innovative medical device utilizes optimum body positioning to achieve the ultimate in local muscle relaxation. This design also provides temporary relief from minor muscular aches due to overexertion.

Dr. Inoue spent decades researching the relationship between body movement and minor muscular pain & tension caused by fatigue or overexertion. He was inspired by watching goldfish swim. Dr. Inoue thought that perhaps the steady undulation causes the fish's whole body to be massaged and provides it with some of its benefits. He wondered what would be the effect on people if they could move like fish?

Engineers took his basic concept and created the technical design for The Chi Machine®. "The pattern of motion that the machine generates in the human body uses minimum effort and achieves maximum results," says Dr. Inoue. "It relaxes muscles locally, temporarily increases our comfort, adds temporary relief to minor tension in muscles from fatigue or overexertion. It can also increase our sense of comfort, through the temporary relief of minor muscle pain and aches, due to fatigue or overexertion."

Regulated as a
Class 1 Medical Device (USA)
Class 2 Medical Device (Canada)

The Chi Machine® vs. an imitation product's
method of creating the swinging motion

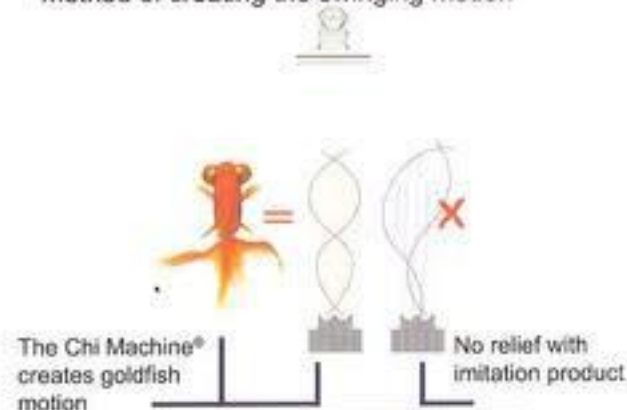

Features & Benefits

Massage is the manipulation of soft tissue for the purpose of increasing overall comfort by temporarily relieving minor muscular pain and discomfort caused by fatigue or overexertion. The first records of massage date back 3000 years to early Chinese folk medicine. Western civilizations were introduced to the benefits of massage by Greek and Roman physicians.

The Chi Machine® blends the ancient arts with modern technology for a unique experience. The Chi Machine® maximizes the movement of the body which helps to provide prompt temporary relief. Fatigue or excess physical exertion can generally lead to minor muscular aches, pains and tension. The Chi Machine® supports temporary relief from these conditions without creating undue stress to bodily structures or additional muscular overexertion.

Because the massage with The Chi Machine® is performed while you are reclining and relaxed, with weight removed from the spine, you will immediately notice a sense of extreme comfort and pleasure.

Since 1990, over Five million Chi Machines® have been sold worldwide. The Chi Machine® has gained a tremendous reputation for its solid construction and outstanding benefits. Unlike many types of medical devices, The Chi Machine® is easy to use, providing excellent results without stress or injury to the body.

How do I use my Chi Machine?

Place the machine on a rug. Lie on your back and place your ankles in the grooves of the footrest on top of the machine. Set the timer switch to activate the machine for the number of minutes you desire. First time users should start with 2-5 minutes.

The footrest will begin to move back and forth, carrying your feet briskly from side to side. This movement will carry up through your legs and torso, throughout your entire body, causing you to wiggle back and forth - that's right, like a fish swimming in water!

All you have to do now is relax and breathe deeply, letting The Chi Machine® rock you gently but firmly for a maximum of 15 minutes.

When the time is up and the machine stops, remain lying down for another two to five minutes, relaxing and breathing deeply, enjoying the remarkable comfort and the temporary feelings of relief.

Although The Chi Machine® is an effective and safe method of massage therapy, its use is NOT advised in the following situations without prior medical advice: during pregnancy, within three months of having had major surgery or bone fracture, if suffering from serious heart disease, a serious infection, any bleeding injury, or if suffering from epilepsy. Additionally, the Therapeutic Massager should not be used within 30 minutes after eating a meal. If severe pain occurs, it should be thoroughly investigated.

The Mid Position



The Goldfish Position



The Stretched Back Position



www.energywellnessproducts.com

**More features of the new SDM-888
The Chi Machine® include:**

- 1) Improved noise insulation - A noise absorbing cushion is added inside to effectively dampen the sound level. Additional stability has been added to the footrest to increase durability and sturdiness.
- 2) With the upgraded material used, The Chi Machine SDM-888 has a smoother and quieter movement.
- 3) A stronger base and construction allows people of all sizes to enjoy the pleasure of The Chi Machine.
- 4) The new fully enclosed design of the motor provides better protection during shipping and transportation. The electrically insulated design of the entire machine ensures the highest safety standard.

**ASK HOW TO EARN YOUR
PRODUCTS FREE!**

**THE BENEFITS OF A
THERAPEUTIC MASSAGER**

- Temporarily relieves minor muscle aches, pains and tension caused by fatigue or overexertion
- Temporarily increases local blood circulation
- Locally relaxes muscles

CALL YOUR REPRESENTATIVE
NOW AND LEARN ABOUT
THE PRODUCTS AND
OPPORTUNITIES
THAT AWAIT YOU!

Try Before You Buy!
• 14 Day Money Back Guarantee •

If you would like to improve your energy level, circulation, oxygenation, lymph drainage, and overall well being, than the Original Chi Machine is for you. You can now get the one that has been copied, but never duplicated in a new and improved version. Some would think it impossible to improve our already great product but we did. With a quieter, smoother, and even more relaxing Chi Machine available, now is the perfect time to invest in your health. Don't wait another day to try our product. We take it with us whenever we travel in the durable custom made carry bag.

John Oldham
President Level Distributor

Working in the wellness field I am continually approached with people wanting me to look at their products. I don't need to because I have been sharing the Chi Machine with others for many years now and I feel that it is the best wellness device on the market. Everyone who uses the Chi Machine gets results. Whether it is pain relief or feeling of well-being, I have seen many positive results. The thing I love most about the Chi Machine is that you don't have to have a medical condition to benefit from it. Taking time every day to relax your mind and your body to the movement of the Chi Machine, focusing on taking deep breathes is a stress reducer and a benefit that definitely leads to better health!

Mara Gerko, CNHP
President Level Distributor

**The
CHI MACHINE®
SDM-888**

"A product that blends the ancient arts with modern technology for a unique experience"

**14 Day
Money Back
Guarantee**

**2
Year
Warranty**



Call today for a free demonstration:

 **HTE USA™**
HTE CANADA™
© 2012 HSIN TEN ENTERPRISES USA, INC. All rights reserved. 0112



www.energywellnessproducts.com